













# Hospital masks



Non Invasive Ventilation (NIV) dedicated masks from emergency to intensive and subacute care

**PHILIPS**

# Patient interface selection

EE versions (transparent elbow with entrainment valve) of all masks below to be used with Single limb circuit with leak connector	Patient characteristics					Clinical status				
	With claustrophobia or anxiety	With facial abnormalities	With no teeth	Easier mouth access	Mouth breathing	Immediate need of ventilation	Acute and severe	Stable and improving	Chronic and stable	Post-surgical / Recovery
<b>Full face mask</b>   Performax*	✓	✓	✓	✓	✓	✓	✓	✓	✓	
<b>Oro-nasal mask</b>   AF811*  PerformaTrak Capstrap*  PerformaTrak*  Image 3*			✓	✓	✓		✓	✓	✓	✓
			✓	✓	✓		✓	✓	✓	✓
			✓		✓		✓	✓	✓	✓
			✓		✓			✓	✓	✓
<b>Nasal mask</b>   PerformaTrak Nasal Capstrap  Contour Deluxe			✓	✓				✓	✓	
			✓	✓				✓	✓	

\* SE version (with blue elbow) also available to be used:

- with dual limb circuit 
- or
- with Single limb circuit with expiratory valve 



**Performax EE**



**Performax SE**

**Performax is intended for patients in acute distress or in need of NIV therapy in intensive and subacute care units**

- Unique contour face design
- Soft and large silicon cushion
- Intuitive and adjustable pre-shaped headgear
- Quick clip snaps
- Compatible with single limb and dual limb circuits

**Part numbers Performax**

PN	Description	Size	Qty / Box
1052554	Performax EE and SE*	L	1 mask and 5 headgear
1052569	Performax SE	L	1 mask and 5 headgear
1052553	Performax EE and SE*	S	1 mask and 5 headgear
1052568	Performax SE	S	1 mask and 5 headgear



**Part numbers Accessories**

PN	Description	Size	Qty / Box
1051816	Performax single use headgear	L	5
1051857	Performax single use headgear	S	5
1045313	Performax exhalation elbow (EE) replacement		1
1043210	Performax standard elbow (SE) replacement		1

\* EE is entrainment valve elbow.  
SE is standard elbow.



**AF811 CapStrap EE**



**AF811 CapStrap SE**

## **AF811 CapStrap the first oro-nasal hospital Gel mask improving comfort for long-term NIV therapy**

- Unique gel cushion with pliable silicone seal which adapts to all morphologies
- Foam forehead pads provide additional comfort and stability
- Capstrap headgear system simplifies mask application, facilitates oral access

### **Part numbers AF811 CapStrap**

PN	Description	Size	Qty / Box
1063015	AF811 Oro-nasal mask with CapStrap EE	L	1
1063016	AF811 Oro-nasal mask with CapStrap EE	M	1
1063047	AF811 Oro-nasal mask with CapStrap EE	S	1
1063051	AF811 Oro-nasal mask with CapStrap EE	L	5
1063052	AF811 Oro-nasal mask with CapStrap EE	M	5
1063053	AF811 Oro-nasal mask with CapStrap EE	S	5
1063048	AF811 Oro-nasal mask with CapStrap SE	L	1
1063049	AF811 Oro-nasal mask with CapStrap SE	M	1
1063050	AF811 Oro-nasal mask with CapStrap SE	S	1
1063054	AF811 Oro-nasal mask with CapStrap SE	L	5
1063055	AF811 Oro-nasal mask with CapStrap SE	M	5
1063056	AF811 Oro-nasal mask with CapStrap SE	S	5
1019543	CapStrap replacement headgear		5
452040	Nasogastric tube sealing pad		10



**PerformaTrak**



**PerformaTrak SE**

**PerformaTrak and PerformaTrak SE are the first intention masks clinically used for patients in intensive or subacute care units**

- Double cushion technology
- Adjustable forehead
- Ball clip headgear system

**Part numbers PerformaTrak**

PN	Description	Size	Qty / Box
1012635	PerformaTrak	L	1
1018591	PerformaTrak	L	10
1012574	PerformaTrak	M	1
1018590	PerformaTrak	M	10
1012573	PerformaTrak	S	1
1018589	PerformaTrak	S	10
1015788	Single use headgear with clips		10

**Part numbers PerformaTrak SE**

PN	Description	Size	Qty / Box
1012641	PerformaTrak SE	L	1
1018595	PerformaTrak SE	L	10
1012640	PerformaTrak SE	M	1
1018606	PerformaTrak SE	M	10
1012639	PerformaTrak SE	S	1
1018605	PerformaTrak SE	S	10
1015788	Single use headgear with clips		10





**PerformaTrak  
Capstrap**



**PerformaTrak  
Capstrap SE**

**PerformaTrak Capstrap and PerformaTrak Capstrap SE are first intention masks for intensive or subacute care**  
**The new advanced headgear system adds stability and comfort compared to the standard headgear**

The new Capstrap system provides the following benefits:



- Enhanced stability and comfort
- Simple application and reapplication
- Easy oral access

#### Part numbers PerformaTrak Capstrap

PN	Description	Size	Qty / Box
1019507	PerformaTrak Capstrap	L	10
1019498	PerformaTrak Capstrap	M	10
1019492	PerformaTrak Capstrap	S	10
1063835	PerformaTrak Capstrap, Single	L	1
1063834	PerformaTrak Capstrap, Single	M	1
1063833	PerformaTrak Capstrap, Single	S	1

#### Part numbers PerformaTrak Capstrap SE

PN	Description	Size	Qty / Box
1019542	PerformaTrak Capstrap SE	L	10
1019541	PerformaTrak Capstrap SE	M	10
1019540	PerformaTrak Capstrap SE	S	10
1063898	PerformaTrak Capstrap SE, Single	L	1
1063897	PerformaTrak Capstrap SE, Single	M	1
1063836	PerformaTrak Capstrap SE, Single	S	1

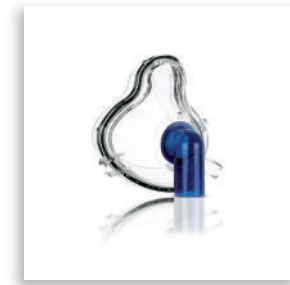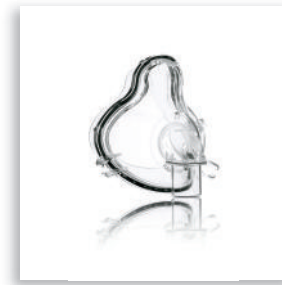

**Image 3**

**Image 3 SE**

**Image 3 and Image 3 SE are the first masks entirely designed for NIV in hospital, for intensive and subacute care units**

- Single cushion technology
- Ball clip headgear system
- Disposable mask

**Part numbers Image 3**

PN	Description	Size	Qty / Box
1018584	Image 3	L	10
1008897	Image 3	L	1
1018583	Image 3	M	10
1004884	Image 3	M	1
1018581	Image 3	S	10
1004879	Image 3	S	1
1004877	Single use 3 straps headgear		1

**Part numbers Image 3 SE**

PN	Description	Size	Qty / Box
1018577	Image 3 SE	L	10
1010879	Image 3 SE	L	1
1018576	Image 3 SE	M	10
1009593	Image 3 SE	M	1
1018575	Image 3 SE	S	10
1009591	Image 3 SE	S	1
1004877	Single use 3 straps headgear		1



**PerformaTrak  
Nasal**



**PerformaTrak  
Nasal Capstrap**

**PerformaTrak Nasal mask is designed for low acuity and clinically stable patients in subacute care units and for post-surgical patients in need of NIV support**

The innovative Capstrap headgear facilitates fitting and improves stability on the patient face, reducing pressure on sensitive facial points.

- Dual flap silicone cushion
- Adjustable forehead
- Capstrap headgear system

**Part numbers PerformaTrak Nasal**

PN	Description	Size	Qty / Box
1048710	PerformaTrak Nasal Capstrap	M/L	10
1048708	PerformaTrak Nasal Capstrap	S	10
1048709	PerformaTrak Nasal Capstrap	P	10





**Contour Deluxe**

**Contour Deluxe is suggested for low acute patients, not mouth-breathing, in subacute care units**

- Dual flap silicone cushion
- Lightweight mask
- Improved headstrap characteristics

**Part numbers Contour Deluxe**

PN	Description	Size	Qty / Box
1016695	Contour Deluxe	M/L	5
1016696	Contour Deluxe	M/L	10
1016697	Contour Deluxe	S	5
1016698	Contour Deluxe	S	10
1016981	Contour Deluxe	P	5
1016983	Contour Deluxe	P	10
1007365	Single use headgear	S/M/L	6
1018472	Single use headgear	P	6



## Part numbers circuits and accessories

PN	Description	Qty / Box
<b>Dry circuits</b>		
582073	NIV circuit BiPAP ventilator without bacterial filter	10
9624C	NIV circuit Focus with bacterial filter and O <sub>2</sub> entrainment port	10
9624	NIV circuit Focus with bacterial filter	10
1065830 <sup>(1)</sup>	NIV circuit with main flow filter and DEP with filter exhalation port and filter	10
1065832 <sup>(1)</sup>	NIV circuit with main flow filter and DEP with filter exhalation port (no filter)	10


PN	Description	Qty / Box
<b>Valves and Filters</b>		
302312	Plateaux exhalation valve	1
332113	Whisper Swivel II exhalation valve	1
1065775 <sup>(2)</sup>	DEP with filter exhalation port (no filter)	10
312149	Bacteria/virus Filter (for use with NIV circuit with DEP filtered exhalation port)	10
1014047	Bacteria/virus Filter (for use with NIV circuit with DEP filtered exhalation port)	1
312149	Exhalation port	10
312010	O <sub>2</sub> enrichment port	10
342077	Bacterial filter low resistance for BiPAP ventilator	10
452040	Naso-gastric tube sealing pad	10

### (1) QTY = 10

Noninvasive patient circuit with main flow bacteria filter, single-use

1065830 With exhalation port filter


1065832 Without exhalation port filter

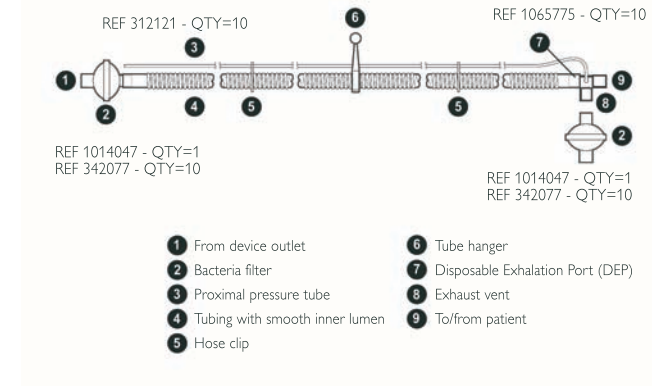

### (2) QTY = 10

Disposable Exhalation Port (DEP)

1065775

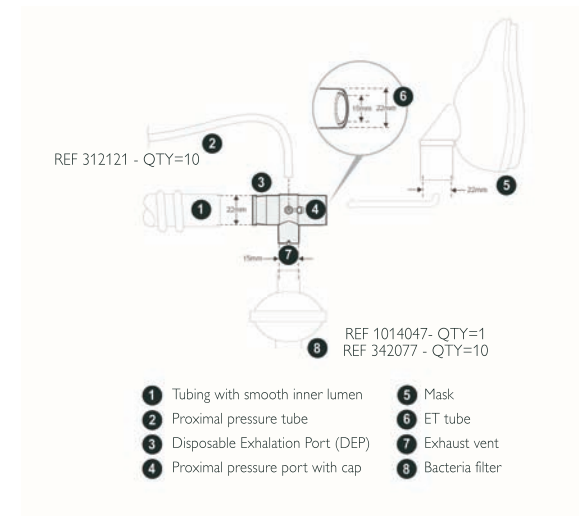





**Philips Healthcare is part of  
Royal Philips Electronics**

**How to reach us**

[www.philips.com/healthcare](http://www.philips.com/healthcare)

[healthcare@philips.com](mailto:healthcare@philips.com)

fax: +31 40 27 64 887

**Asia**

+852 2821 5888

**Europe, Middle East, Africa**

+49 7031 463 2254

**Latin America**

+55 11 2125 0744

**North America**

+1 425 487 7000

800 285 5585 (toll free, US only)

**Customer Service**

+1 724 387 4000

800 345 6443 (toll free, US only)

**Respironics Europe, Africa, Middle East**

+33 1 47 52 30 00

**Respironics Germany**

+49 8152 93 06 0

**Respironics France**

+33 2 51 89 36 00

**Respironics Italy SRL**

+39 039 203 6667

**Respironics UK**

+44 800 1300 845

**Respironics Switzerland**

+41 62 745 1750

**Respironics Asia Pacific**

+852 3194 2280

**Respironics Australia**

+61 2 9666 4444

[www.philips.com/respironics](http://www.philips.com/respironics)



© 2009 Koninklijke Philips Electronics N.V.

All rights are reserved.

Respironics, Performax, PerformaTrak, Capstrap, Total Face; Contour Deluxe, BiPap, Vision and Focus are trademarks of Respironics, Inc. and its affiliates.

Philips Healthcare reserves the right to make changes in specifications and/or to discontinue any product at any time without notice or obligation and will not be liable for any consequences resulting from the use of this publication.

Printed in France.

RIT HOSP 029 E \* DEC 2009