


ASAP^{LP}TM

ASAP[®]

ASPIRATION CATHETERS AND KITS


Your Options for
Effective and Quick
Clot Aspiration

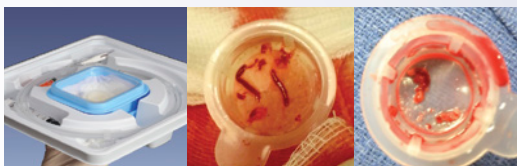

ASAPLP™ AND ASAP® ASPIRATION CATHETERS KEY FEATURES



ASAPLP™
ASPIRATION CATHETER KIT

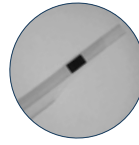

EFFECTIVE

- 100% stainless steel braided, kink resistant 145 cm low profile aspiration catheter
- Improved tip radiopacity due to radiopaque platinum iridium marker band
- Bright colored wire exchange lumen

QUICK ASPIRATION

- Fits into 6F guide catheter of ≥ 0.066 " with ease
- A 4F catheter designed for maximum aspiration
- Long 30 cm hydrophilic coating promoting deliverability
- Longer 20 cm rapid exchange lumen designed for better wire trackability

OPTIONS

- Packaged with stiffening stylet



ASAP®



EFFECTIVE

- Single extrusion, 100% stainless steel braided, kink resistant 140 cm aspiration catheter
- Tapered tip for smooth catheter transition

QUICK ASPIRATION

- Large aspiration lumen to facilitate quick aspiration
- 3 depth-positioning nonradiopaque markers at 90 cm, 100 cm, 110 cm
- All-in-one kit with all the necessary components to perform a catheter aspiration procedure

OPTIONS

- Standalone MiniStop® waste basin with lid and 70-micron filter



CONVENIENCE AT YOUR FINGERTIPS


Both catheters are available in the ASAP Kit with all the necessary components to perform a catheter aspiration procedure safely, quickly, and effectively. With our all-in-one tray, simply remove the lid to begin the procedure.

KIT COMPONENTS


- (1) Rapid exchange kink resistant ASAPLP or ASAP aspiration catheter
- (2) Merit VacLok® 30 mL aspiration syringes
- (1) 4 mL RXP® flushing syringe
- (2) Filter baskets, 70 micron
- (1) MicroStop™ closed waste basin, splatter proof, 60 mL fluid capacity
- (1) 21.5 cm tubing set with stopcock

A FAST, HANDY SOLUTION

The ASAP kits are designed for the quick removal of fresh, soft emboli and thrombi from vessels of the arterial system.

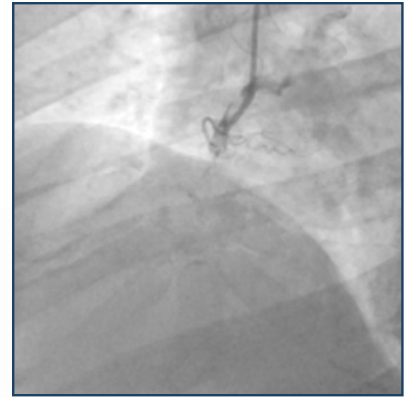

Pre Aspiration


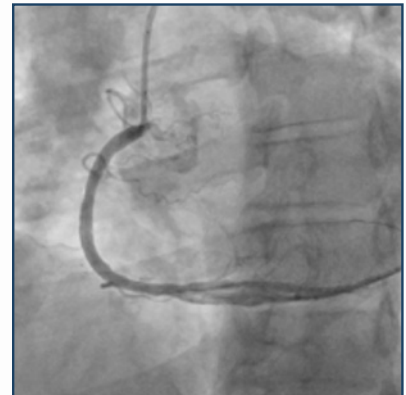


Post Aspiration





ASAP[®]

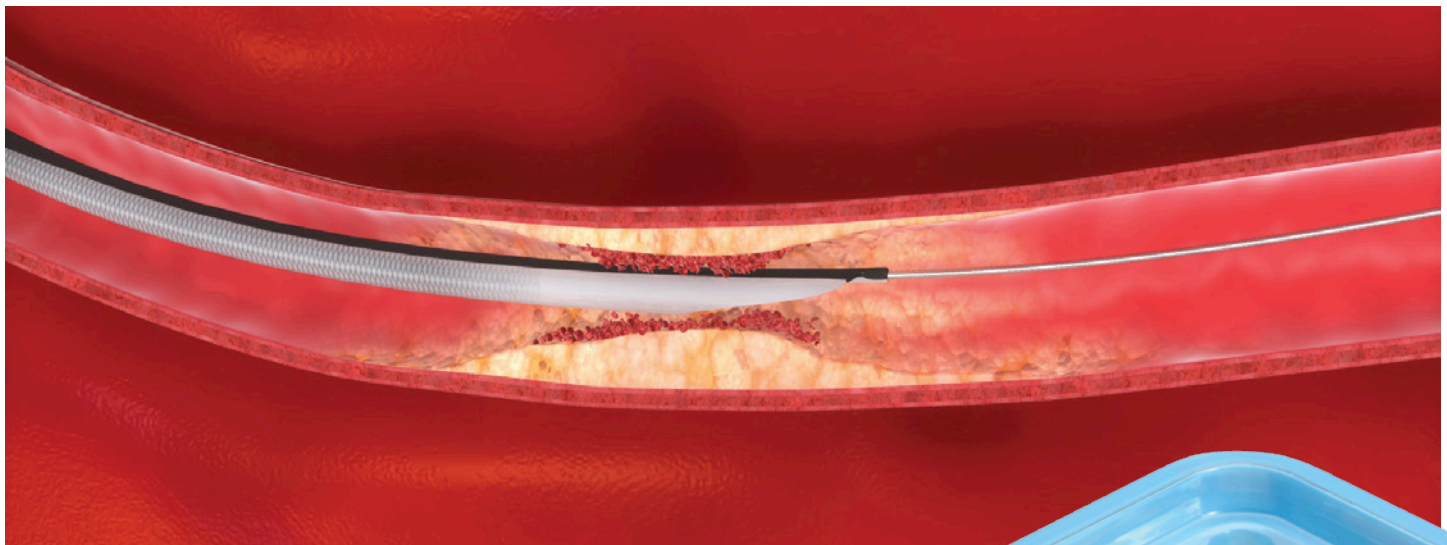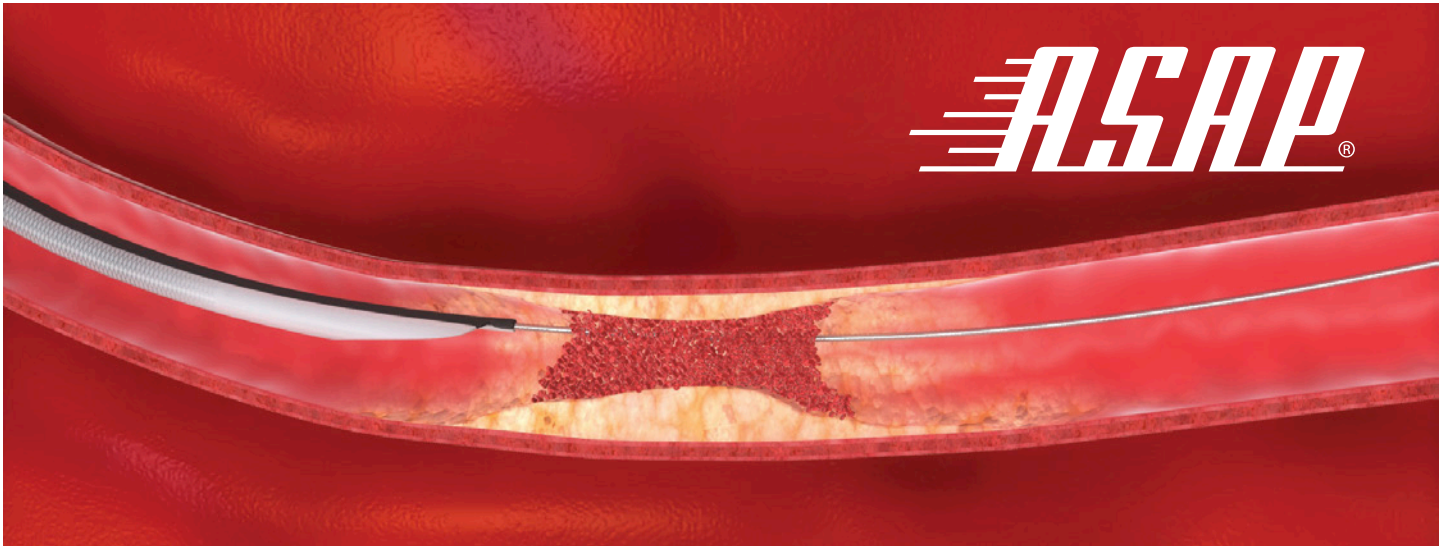

ORDERING INFORMATION

Catalog Number	Description	Box Quantity
ASAP100	ASAP Aspiration Catheter & Kit	1
ASAPLP	ASAPLP Catheter & Kit	1
MMSF200	MiniStop Waste Basin with Lid and Filter	20

U.S. Patent Nos. 6,077,698; 6,278,018; 6,603,040; 6,924,390; 7,138,541



Merit Medical Systems, Inc. • 1600 West Merit Parkway • South Jordan, Utah 84095 • 1-801-253-1600 • 1-800-35-MERIT

Merit Medical Europe, Middle East, and Africa (EMEA) • Amerikalaan 42, 6199 AE Maastricht-Airport • The Netherlands • Tel: +31 43 358 82 22

Merit Medical Ireland Ltd. • Parkmore Business Park West • Galway, Ireland • +353 (0) 91 703 733

Free phone for specific country: Austria 0800 295 374 • Belgium 0800 72 906 (Dutch) 0800 73 172 (French) • Denmark 80 88 00 24 • France 0800 91 60 30

Finland 0800 770 586 • Germany 0800 182 0871 • Ireland (Republic) 1800 553 163 • Italy 800 897 005 • Luxembourg 8002 25 22 • Netherlands 0800 022 81 84

Norway 800 11629 • Sweden 020 792 445 • UK 0800 973 115

www.merit.com

401951010/A ID 053013